


Skills Checklist : Measuring Temperature

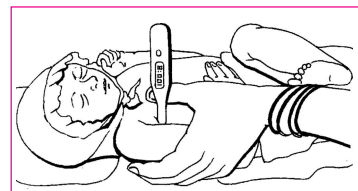

Skills Checklist	For Peer Record					For Trainer Record					
	1	2	3	4	5	6	7	8	9	10	11
1) Take thermometer out of its storage case, hold at broad end, and clean the shining tip with cottonball soaked in spirit.											
2) Press the pink button once to turn the thermometer on. You will see "188.8" flash in the center of the display window, then a dash (-) then the last temperature taken, and then three dashes (- - -) and a flashing "F" in the upper right corner.											
3) Hold the thermometer upward and place the shining tip in the center of the armpit. Place arm against it. Do not change the position.											
4) You will hear a beep sound every 4 seconds while the thermometer is recording the temperature. When you hear 3 short beeps, look at the display. If "F" is flashing, and numbers are changing, leave the thermometer in place. When "F" stops flashing and the numbers stop changing, remove the thermometer.											
5) Read the number in the display window.											
6) Record the temperature reading on the form											
7) Turn the thermometer off by pushing the pink button one time.											
8) Clean the shining tip of the thermometer with a cotton ball soaked in spirit.											
9) Place thermometer back in its storage case.											

Name of the ASHA : _____

Date : _____

Name of the trainer : _____

Total score : _____

Block : _____